

## Bio-Adembeads Antibodies Goat anti-Mouse IgG G0503

For research use only

### INTRODUCTION

Bio-Adembeads Antibodies Goat anti-Mouse IgG were developed for the positive selection or depletion of cells labelled with primary Mouse IgG antibodies. They can also be utilized for the positive selection or depletion of subcellular material, bacteria or other microorganisms labelled with primary Mouse IgG antibodies. They can also be used for the positive selection or depletion of Mouse B cells expressing surface IgG.

### PRODUCT DESCRIPTION

Bio-Adembeads Antibodies Goat anti-Mouse IgG are monodispersed superparamagnetic particles coated with affinity purified Goat anti-Mouse IgG covalently bound to the surface. These secondary antibodies are Fc<sub>γ</sub> Fragment Specific and bind all mouse IgG subclasses. Cross reactivity to Human, bovine and Horse serum proteins is minimal. They are produced under aseptic conditions and are supplied in an aqueous suspension containing 0.05% Proclin 300.

### PHYSICAL CHARACTERISTICS

Diameter: 300 nm (CV max 20%)

Magnetic susceptibility: approx. 40 emu/g

Specific surface area: 10 m<sup>2</sup>/g

Iron oxide content: approx. 70%


### PRINCIPLE OF BIO ADEMBEADS ANTIBODIES SEPARATION

#### Indirect Coating

Cells are labeled with primary Mouse IgG antibody. Subsequently, cells are rosetted with Bio-Adembeads Antibodies Goat anti-Mouse IgG for a magnetically labeling.

#### Direct Coating

In this case the primary Mouse antibody of the IgG-class is bound firstly to the Bio-Adembeads Antibodies Goat anti-Mouse IgG before incubation with the cell suspension.

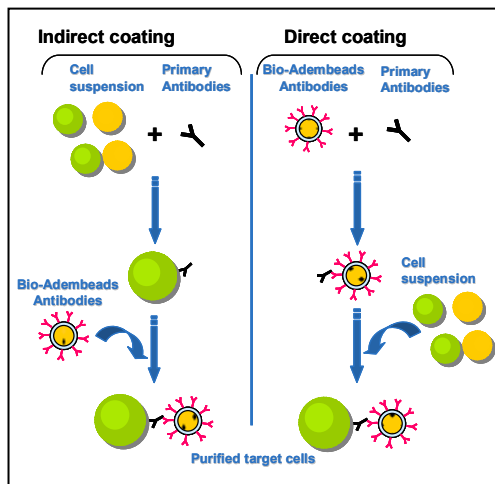

#### Note:

**We recommended that you use indirect technique for cell separation and other formats of purification.**

**The type and amount of antibody required for optimal coating and performance of the beads will vary with the mAb-affinity and the antigen density on the surface. The optimal coating of beads to achieve an efficient binding of cells should be determined in each individual case by titration.**

### INSTRUCTIONS FOR USE

#### A) Indirect coating for cell separation

Before cells can be isolated with Bio-Adembeads Antibodies Goat anti-Mouse IgG in the indirect technique, the cells must be incubated with specific Mouse primary antibodies against surface markers on the target cells.

1. Incubate the heterogeneous cell suspension with the primary Mouse IgG antibodies for 30 minutes at 2-8°C providing gentle mixing. Each primary antibody should be added in excess to occupy all antigen binding sites on the cell surface.
2. Collect the incubated cells by centrifugation at 800g for 10 minutes. Discard supernatant.
3. Resuspend and wash the incubated cells 2 times with Hanks Balanced Salt Solution, pH 7.4 to remove all unbound antibody. The cells are now covered with specific primary Mouse antibodies. A cell concentration of 10-50 mill/ml sample is recommended.
4. Add minimum 30 µl of Bio-Adembeads Antibodies Goat anti-Mouse IgG per ml of primary antibody covered cells. (Titration ensures an optimal use of the Bio-Adembeads products)
5. Incubate for approx 15-30 minutes with titling and rotation mixer at 2-8°C.
6. Place the tube in a Adem-Mag SV or MV for 2-5 minutes and pipette off the supernatant. Alternatively, incubation cells medium can be added after the incubation to the height of the magnet to avoid trapping of non-bound cells.
7. Positively selected cells attached to the Bio-Adembeads should be washed four to five times with cells culture medium.

#### B) Direct coating for cell separation

1. Resuspend the Bio-Adembeads Antibodies Goat anti-Mouse IgG by pipetting and vortexing. Avoid foaming.
2. Pipette the volume to be used into the desired test tube.
3. Place the tube in a magnet (see Related Product) for 1min.
4. Pipette off the supernatant carefully, leaving beads undisturbed.
5. Remove the test tube from the magnet (see Related Product) and resuspend the beads

carefully in the original sample volume with adequate buffer (Immobilisation Buffer recommended, see Related Product).

6. Resuspend thoroughly by vortexing and/or shaking.
7. Add maximum 1-5 µg Mouse IgG to 100 µl of washed Bio-Adembeads Antibodies Goat anti-Mouse IgG.
8. Incubate the mixture from 5 to 30 minutes at 20°C (incubation temperature depends on antibody stability) with titling and rotation for even mixing.
9. Place the tube in the magnet (Adem-Mag SV or MV, see related product) and leave to separate for 2 minutes.
10. Discard the supernatant while the tube is on the Adem-Mag magnet, taking care not to disturb the Adembeads. Remove the tube from the Adem-Mag magnet and resuspend the beads in the adequate buffer (Immobilisation Buffer). Mix gently for a few seconds, before applying the magnet again.
11. Repeat step 9-10 twice but the second time resuspend the beads in Hanks Balanced Salt Solution or incubation cells medium.
12. Incubate precoated Bio-Adembeads Antibodies Goat anti-Mouse IgG and the heterogeneous cell suspension for approx. 15-30 minutes with titling and rotation at 2-8°C. Use a final Bio-Adembeads concentration suitable for your experiment.
13. Place the tube in a Adem-Mag SV or MV for 2-5 minutes and pipette off the supernatant. Alternatively, incubation cells medium can be added after the incubation to the height of the magnet to avoid trapping of non-bound cells.
14. Positively selected cells attached to the Bio-Adembeads should be washed four to five times with cells culture medium.

## ADDITIONAL MATERIAL REQUIRED

- Magnetic device
- Rotation device
- Test tubes
- Related products: Buffer solutions
  - *Immobilisation Buffer (# 10301)*Magnetic Devices
  - Adem-Mag SV, 1.5 ml (# 20101)
  - Adem-Mag MV, 15 ml (# 20102)
  - Adem-Mag HV, 50 ml (# 20103)

## STORAGE/STABILITY

When stored in unopened vials at 2-8°C, Bio-Adembeads are stable until expiration date printed on the label.

The Bio-Adembeads must be maintained in liquid during storage and all handling steps. Drying will result in reduced performance. Do not freeze the product.

## PRECAUTIONS

Precautions should be taken to prevent bacterial contamination of protein-coated Adembeads. If cytotoxic preservatives are added they must be carefully removed before use by washing.

## WARNINGS AND LIMITATIONS

**For in vitro research only.** Not for use in human diagnostic or therapeutic procedures.

Proclin 300 is toxic if ingested. **Avoid pipetting by mouth.**

## WARRANTY

The products are warranted to the original purchaser only to conform to the quality and contents stated on the vial and outer labels for duration of the stated shelf life.

Ademtech's obligation and the purchaser's exclusive remedy under this warranty is limited either to replacement, at Ademtech's expense, of any products which shall be defective in manufacture, and which shall be returned to Ademtech, transportation prepaid, or at Ademtech's option, refund of the purchase price.

Claims for merchandise damaged in transit must be submitted to the carrier.