



Smart D-N-Adem-Kit

For blood and PBMC

(cat #06111)

**Instruction manual for
gDNA capture and direct amplification**

■ **Blood**

■ **PBMC**

Ademtech SA
Parc scientifique Unitec 1
4, allée du Doyen G. Brus
33600 PESSAC
France

Tel: +33557020201
Fax +3355702020

Visit our Web site: www.ademtech.com

Table of contents

General Overview 4

Smart D-N-Adem-Kit for Blood and PBMC 4

Product Description 4

Smart D-N-Adem-Kit Protocol 6

- Smart D-N-Adem-Kit procedure 6

- General Guidelines 6

- Protocol for gDNA isolation from blood and PBMC 7

Troubleshooting 11

Warranty 11

General Overview

Biomagnetic separation technology is a simple technique based on the separation of **superparamagnetic beads** using a magnetic field. When added to a complex medium, the magnetic particles will bind to the target. This interaction is based on the specific affinity of the ligand to the surface of the beads. The resulting target-bead complex can be removed from the suspension using a magnet. The inherent benefits of magnetic handling allow for easy washing, separation and concentration of the target without any need of centrifugation or columns.

Superparamagnetic beads exhibit magnetic properties only when placed within a magnetic field and show no residual magnetism when removed from this field.

Smart D-N-Adem-Kit for blood and PBMC

1. Product Description

1.1. General Description

The Smart D-N-Adem-Kit is a new system which is specially designed for the capture of gDNA and direct amplification on magnetic beads. Smart-Adembeads have a specific and proprietary polymer on their surface, designed for DNA capture by electrostatic interactions and compatible with direct PCR and Real-Time PCR without the need to perform an elution step.

PBMC or blood are mixed with the Lysis Buffer. Smart-Adembeads bind specifically to the gDNA. Proteins and other contaminants are then eliminated in the washing step. The purified gDNA bound to the Smart-Adembeads can be used directly for PCR and Real-Time PCR analyses. The Smart D-N-Adem-Kit procedure allows cleaning gDNA in less than 10 minutes, lysis step included. DNA isolation is achieved without phenol, ethanol, chloroform and ionic chaotropes; thus the purified gDNA bound to the Smart-Adembeads demonstrates improved downstream performance in PCR and qPCR. Unlike other purification systems, no elution of the DNA is required, making it possible to maximize the amount of templates available for the reaction achieving greater sensitivity. This method makes this Kit ideal for processing small amount of samples where maximum DNA recovery is critical.

1.2. Kit capacities

Sample	Amount of starting material	Number of isolation
PBMC	10- 5x10 ⁵	100 isolations
Blood	0.1-50µl	100 isolations

Table 1: Number of isolation per amount of starting material

1.3. Reagents provided with the kit

The smart D-N-Adem-Kit includes reagents for performing 100 gDNA isolations. The provided reagents are listed below.

	Amount	Component	Storage
R1	0.5ml	Smart-Adembeads	+ 4°C
R2	0.1 ml	RNase A (2mg/ml)	+ 4°C
R3	0.5 ml	Proteinase K (10mg/ml)	+ 4°C
R4	15ml	Lysis Buffer	+ 4°C
R5	15 ml	Binding Buffer	+ 4°C
R6	20ml	Washing Buffer	+ 4°C
R7	5 ml	Amplification Buffer	+ 4°C

Table 2: Reagent provided with the kit

Properly stored Kits are guaranteed until the expiration date. Note that the shipping is realized at room temperature which will not affect its stability. All the components of the kit have been prepared under nucleases free conditions and have been thoroughly tested to ensure optimal performance.

1.4. Required equipment (not supplied as part of the kit)


- Adem-Mag SV (#20101) or MSV (#20104):
- Nuclease free microtubes
- Nuclease free tips
- Disposable gloves



Adem-Mag MSV

2. Smart D-N-Adem-Kit Protocol

2.1. Smart D-N-Adem-Kit procedure

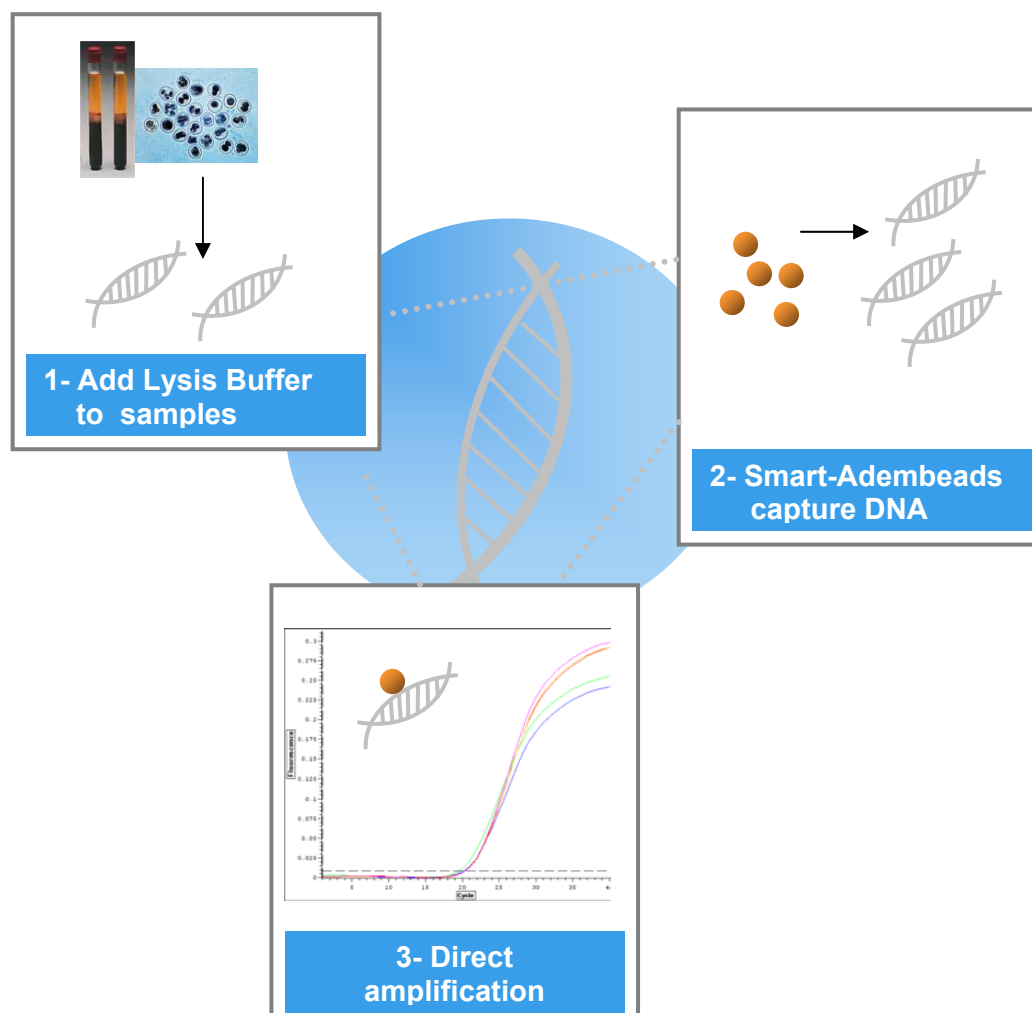

2.2. General Guidelines

- The Smart D-N-Adem-Kit is suitable only for capture and direct amplification of gDNA on magnetic beads from small amounts of samples without the need for an elution step.
- Using UV spectrophotometric measurements (A260/A280 nm), is not recommended as this method is inaccurate for low DNA.
- Before starting the gDNA extraction procedure, carefully read the paragraph **2.3.6**.
- Before starting the gDNA extraction procedure all buffers should be at room temperature to attain optimal performances.

2.3. Protocol for gDNA isolation from blood and PBMC

2.3.1 Sample Preparation:

From whole blood: This procedure is described up to 50µl

- 1- Transfer the whole blood to a 1.5ml microcentrifuge tube.
- 2- Add 50µl or 100µl of Lysis Buffer (see Table 3).
- 3- Add 5µl of Proteinase K and 1µl of RNase. Mix by pipetting (or flick slowly the tube) and incubate at room temperature for 5 minutes (do not exceed 5 minutes).
- 4- Go to paragraph **2.3.2**.

Note: Be careful not to heat the blood.

From PBMC: This procedure is described up to 5×10^5 cells

1. Collect the PBMC and transfer them to a 1.5ml microcentrifuge tube in order to have them in a minimal volume of culture media or PBS (up to 1/5 of the Lysis Buffer volume, see Table 4).
2. Add 50µl to 150µl of Lysis Buffer (see Table 4).
3. Add 5µl of Proteinase K and 1µl of RNase. Mix by pipetting (or flick the tube) and incubate at room temperature for 5 minutes.
4. Go to paragraph **2.3.2**.

Note:

- Refer to Table 3 and Table 4 for the recommended amount of reagents.
- To save time, a solution containing Lysis Buffer, Proteinase K and RNase can be prepared to add directly to the starting samples. Caution, RNase should be added in the last to avoid its early degradation by Proteinase K.

2.3.2 DNA Capture

- 1- Add 1 μ l or 5 μ l of homogenized Smart-Adembeads (see Table 3 or Table 4).
- 2- Add 50 μ l to 150 μ l of Binding Buffer and homogenize by pipetting (see Table 3 or Table 4). Incubate at room temperature for 1 minute.

Note:

- It is important not to increase the incubation time, that could decrease the efficiency of DNA capture.
- To save time, a solution containing the Smart-Adembeads and Binding Buffer should be prepared to add directly to your lysis samples.

2.3.3 Wash [DNA Smart-Adembeads] complexes

- 1- Place the tube on the magnet for at least 1 minute or until supernatant clearing and discard the supernatant. Remove the tube from the magnet and add 100 μ l of Washing Buffer and mix by pipetting.
- 2- Repeat once step 1 using 100 μ l of Washing Buffer.

Note: When removing the supernatant, start by pipetting the potentially formed foam.

2.3.4 Resuspension [DNA Smart-Adembeads] complexes

- 1- Place the tube on the magnet for at least 1 minute or until supernatant clearing and discard slowly the supernatant.
- 2- Resuspend beads in 10 or 50 μ l of Amplification Buffer and mix by pipetting (see Table 3 or Table 4).

Note:

- You may observe the formation of aggregates due to DNA capture by beads, which can be easily resuspended by pipetting.
- If you have drops of blood on the surface of the microtube, transfer the final solution in a new tube.

- 3- Use the final solution directly for PCR and Real-Time PCR (see recommendations in paragraph 2.3.6) or/and store it at -20°C (long-term storage).

2.3.5 Choose your own starting point

Whole Blood	Up to 10µl	10-50µl
Lysis Buffer	50µl	100µl
Proteinase K	5µl	5µl
RNase	1µl	1µl
Smart-Adembeads	5µl	5µl
Binding Buffer	50µl	100µl
Washing Buffer (x2)	100µl	100µl
Amplification Buffer	50µl	

Table 3: Buffer amounts for gDNA preparation

PBMC number	*10 - 10 ²	10 ² – 25x10 ³	25x10 ³ - 10 ⁵	10 ⁵ – 5x10 ⁵
Lysis Buffer	50µl	50µl	100µl	150µl
Proteinase K	5µl	5µl	5µl	5µl
RNase	1µl	1µl	1µl	1µl
Smart-Adembeads	1µl	5µl	5µl	5µl
Binding Buffer	50µl	50µl	100µl	150µl
Washing Buffer (x2)	100µl	100µl	100µl	100µl
Amplification Buffer	10µl	50µl	50µl	50µl

Table 4: Buffer amounts for gDNA preparation from PBMC

*For very small starting sample (less than 10 cells), resuspending the complexes DNA/beads in 5µl of amplification buffer and using all of the solution for the PCR detection is recommended.

2.3.6 Recommendations for direct DNA amplification on Smart-Adembeads

■ General guidelines for PCR

- 1- Homogenize the final solution thoroughly (gDNA/Beads complexes in Amplification Buffer).
- 2- Add the volumes of the previous solution (see below) to the PCR mixtures.
- 3- Perform your PCR reactions following manufacturers directions.
- 4- Make sure the number of cycles is sufficient and the primers have been tested and are satisfactory quality for amplification of small quantities of DNA.

■ Performing Standard PCR:

Using 1-5µl of the final solution (gDNA/Beads complexes in Amplification Buffer) is recommended for 25 or 50µl of PCR reaction. No more than 10µl are suitable.

■ Performing Real-Time PCR for detection:

From whole blood and PBMC, using 1-5µl of the final solution (gDNA/Beads complexes in Amplification Buffer) is recommended for a 25µl of PCR reaction or 1-10µl of the final solution (gDNA/Beads complexes in Amplification Buffer) for a 50µl of PCR reaction.

■ Performing Real-Time PCR for quantification:

From whole blood, the starting quantity of sample must not exceed 0.5µl of whole blood. We recommend diluting your starting sample in PBS (for example : dilute 5µl of whole blood in 45µl of PBS and use 5µl of this solution for the process).

From PBMC, the starting quantity of sample should not exceed 5×10^3 cells. We recommend diluting your starting sample in PBS or culture media. Be careful, the diluted starting sample must not exceed 1/5 of the Lysis buffer volume.

To quantify DNA, the following volumes are recommended,

*using 1-3µl of the final solution (gDNA/Beads complexes in Amplification Buffer) for a 25µl of PCR reaction.

*using 1-5µl of the final solution (gDNA/Beads complexes in Amplification Buffer) for a 50µl of PCR reaction.

Troubleshooting

No or low amounts of PCR products detected after standard PCR or qPCR.

- **DNA degradation:** Ensure that the process is completed in a nuclease free environment to avoid introducing any nucleases during the procedure or later handling.
- **Optimize the lysis step:** From PBMC, incubation time should be increased an additional 5 or 10 minutes.
- **Optimize the washing step:** From whole blood, make sure beads are resuspended correctly in the Washing Buffer in order to eliminate haemoglobin which is a PCR inhibitor.
- **Too much gDNA/beads complexes used as PCR templates:** Decrease the volume of the final solution (gDNA/beads in Amplification Buffer) that is introduced in PCR reactions.
- **Optimize the qPCR or PCR assay:**
 - For results improvement, the first PCR step (initial denaturation) could be increased at 3 to 5 minutes.
 - Optimization parameters include primer design, primer concentration, probe design and probe concentration.
 - For best results, design primers that produce amplicons <150 bp in length.
 - Increase the number of PCR cycles.
 - Determine the optimal primer binding temperature for the PCR cycling program.
 - For long strength fragment, a more processive DNA polymerase is recommended.

Warranty

This product is only for use in research. The purchaser is responsible to validate the performance of this product for any particular use, and to use the product in compliance with any applicable regulations. The products are warranted to the original purchaser only to conform to the quality and contents stated on the vial and outer labels for duration of the stated shelf life. Ademtech's obligation and the purchaser's exclusive remedy under this warranty is limited either to replacement, at Ademtech's expense, of any products which shall be defective in manufacture, and which shall be returned to Ademtech, transportation prepaid, or at Ademtech's option, refund of the purchase price. Claims for merchandise damaged in transit must be submitted to the carrier.